
Supplementary file S1

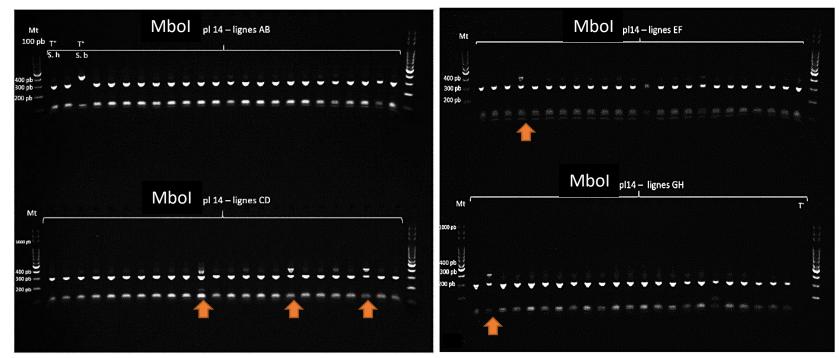

Gel revelation after PCR-RFLP of ITS2 gene. The PCR primers were : Forward: 5'-GGCTGCAGCGTTAACCATTA-3' and reverse: 5'-ACACACCATCGGTACAAA-3'. The amplified fragments were digested with the Mbol enzyme. After digestion the expected fragments are 44, 82, and 379 bp for *S. haematobium*, and fragments of 44, 82, 98 and 281 bp for *S. bovis*. Only bands higher than 100 bp are visible on the gels.